



EASYCAP

EEG Recording Caps and Related Products


EASYCAP GmbH
Am Anger 5
82237 Würthsee
GERMANY

Tel +49 8153 88702-00
www.easycap.de
info@easycap.de

Cap for 96Ch actiCAP slim Standard Layout

Electrode Names and Number Labels

actiCAP
CCCSlim

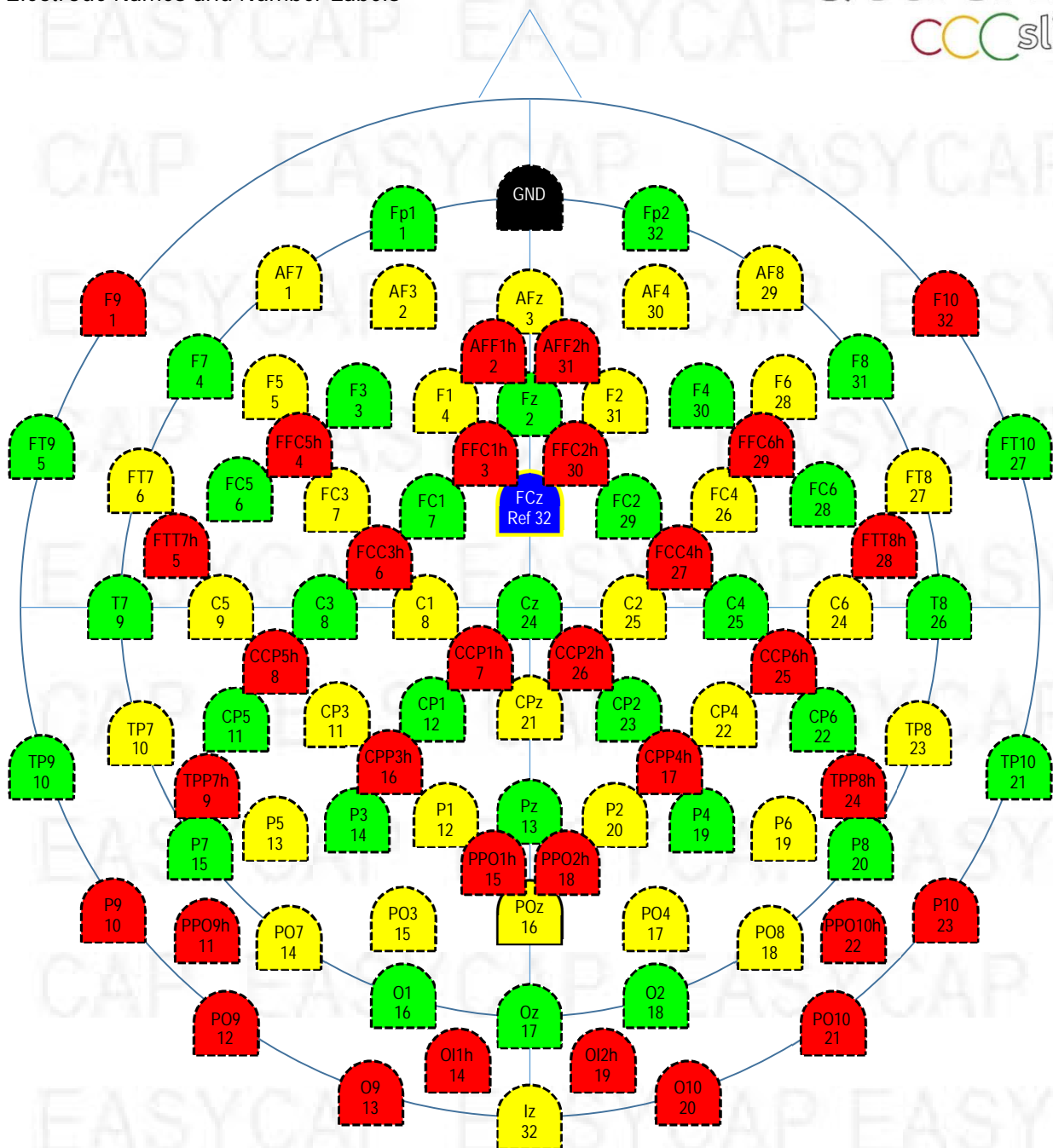

Details

Ordering Information

For ordering please give Article Number, Cap Cut, and Size (e.g. CAS-96, C-Cut, 56):

- Article Number: CAS-96
- Cap Cut: C-Cut or A-Cut
- Size (given in cm head circumference):
Adult caps: 52, 54, 56, 58, 60, 62, 64 (average male: 58, average female: 56)

The catalogue-number comprises the cap with 64 + 2 openings for actiCAP slim electrodes as described. For further information about accessories or consumables, please visit our website or contact our local distributor.

Cap

Standard: White Subinon Cap with integrated chin belt

Sizes 52 – 64 made from High Precision Fabric, Sizes 50 and smaller made from High Comfort Fabric


Options: C-Cut or A-Cut, Size. For further variations, contact us

Insertion of Electrodes

Button the electrodes directly into the cap fabric, as indicated in the figure: Electrodes 1-32, 33-64, and 65-96 are actiCAP-Slim electrodes, terminating at a splitter box for connection to actiCHamp or actiCAP control box. REF and GND (160cm) come with individual connectors.

The velcro straps attached to the cap help you to organize and fixate the cable wires.

For alternative usage (e.g. EOG), electrodes can also be attached to the skin with washers (= double-sided adhesive rings).

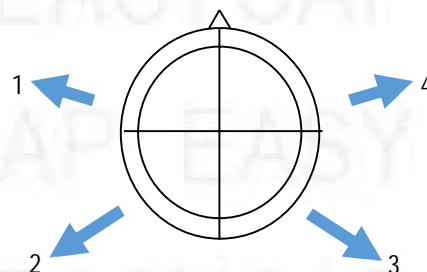

Please note: Each 32Ch electrode set is labelled 1-32; however, the splitter boxes are numbered 1 through 3.

Please use the electrodes of splitter box 1 for channels 1-32 (green positions in the figure).

Please use the electrodes of splitter box 2 for channels 33-64 (yellow positions in the figure).

Please use the electrodes of splitter box 3 for channels 65-96 (red positions in the figure).

If your hardware requires a REF electrode (e.g. actiCAP control box), button REF into the opening at FCz. You may then use Iz position for yellow channel 32.

If your hardware uses NO REF electrode (e.g. actiCHamp), leave the REF electrode unused. You may use FCz, or Iz, for yellow channel 32.

Hints when performing TMS and EEG simultaneously

To minimize the TMS artefacts onto the EEG signal, electrode cables should be led away from the TMS stimulation point. To be able to do so, the electrode cables are not attached to the cap and can be rotated 360° in their cap holes. You may secure a certain cable geometry with the enclosed cable clips and velcro straps. Depending on the site of the TMS stimulation, it may be advisable to move the electrode used as REF away from the stimulation point to another position further away.